

ABRAXIS® Microcystins 0.3 ppb Strip Test Finished Drinking Water 520046

1. Collect Sample




SAMPLE

Collect 1 to 2 mL of sample.


2. Transfer

Squeeze the upper bulb to draw sample into the pipette. The required amount of sample will fit into the stem of the pipette. Any excess will flow into the overflow bulb. Squeeze the upper bulb again to dispense the required amount of sample into conical tube.

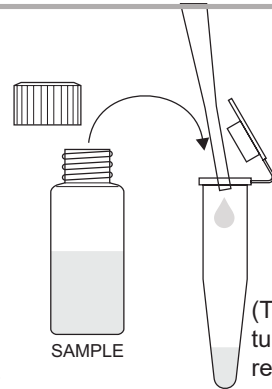

(The conical, flip-top tube contains dried reagents.)

3. Shake and Incubate

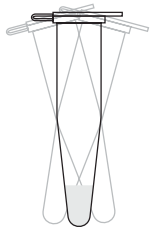

Close the conical, flip-top tube and shake for 30 seconds.

(Dried reagents will dissolve, turning the sample purple.)



Incubate for 10 minutes.



4. Test

Insert test strip into conical, flip-top tube with arrow pointing down.



Incubate for 10 minutes.




5. Dry

Remove test strip. Lay flat and allow to continue developing for 5 minutes.



6. Interpret

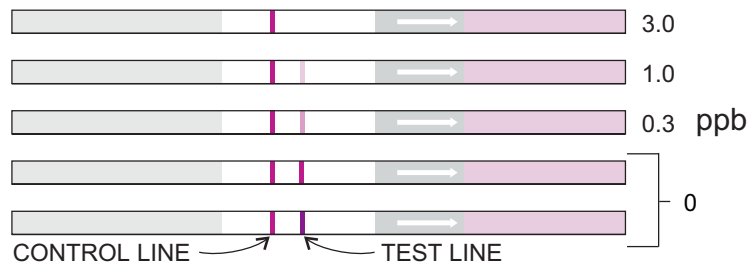

INTERPRET TEST

CONTROL LINE	TEST LINE	INTERPRETATION
NO CONTROL LINE PRESENT	NO TEST LINE PRESENT	INVALID RESULT
CONTROL LINE PRESENT	NO TEST LINE PRESENT	>3 ppb
CONTROL LINE PRESENT	MODERATE INTENSITY TEST LINE PRESENT	BETWEEN 0 AND 3 ppb