

KIT COMPONENTS SUPPLIED BY CUSTOMERS



Tubes for Sampling

****Do NOT re-use any tubes for multiple samples due to a high risk of toxicity and contamination. Dispose of the tubes after single use****


KIT COMPONENTS SUPPLIED BY SENSOREAL



Conjugate Tubes




Diluent Powder

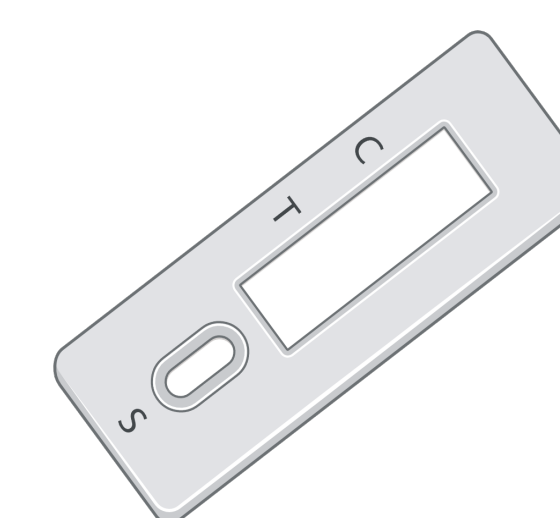


Lateral Flow Test Strip

PART 1: SAMPLE PREPARATION


1. To a tube "D" add **5 mL** of your sample and mix.

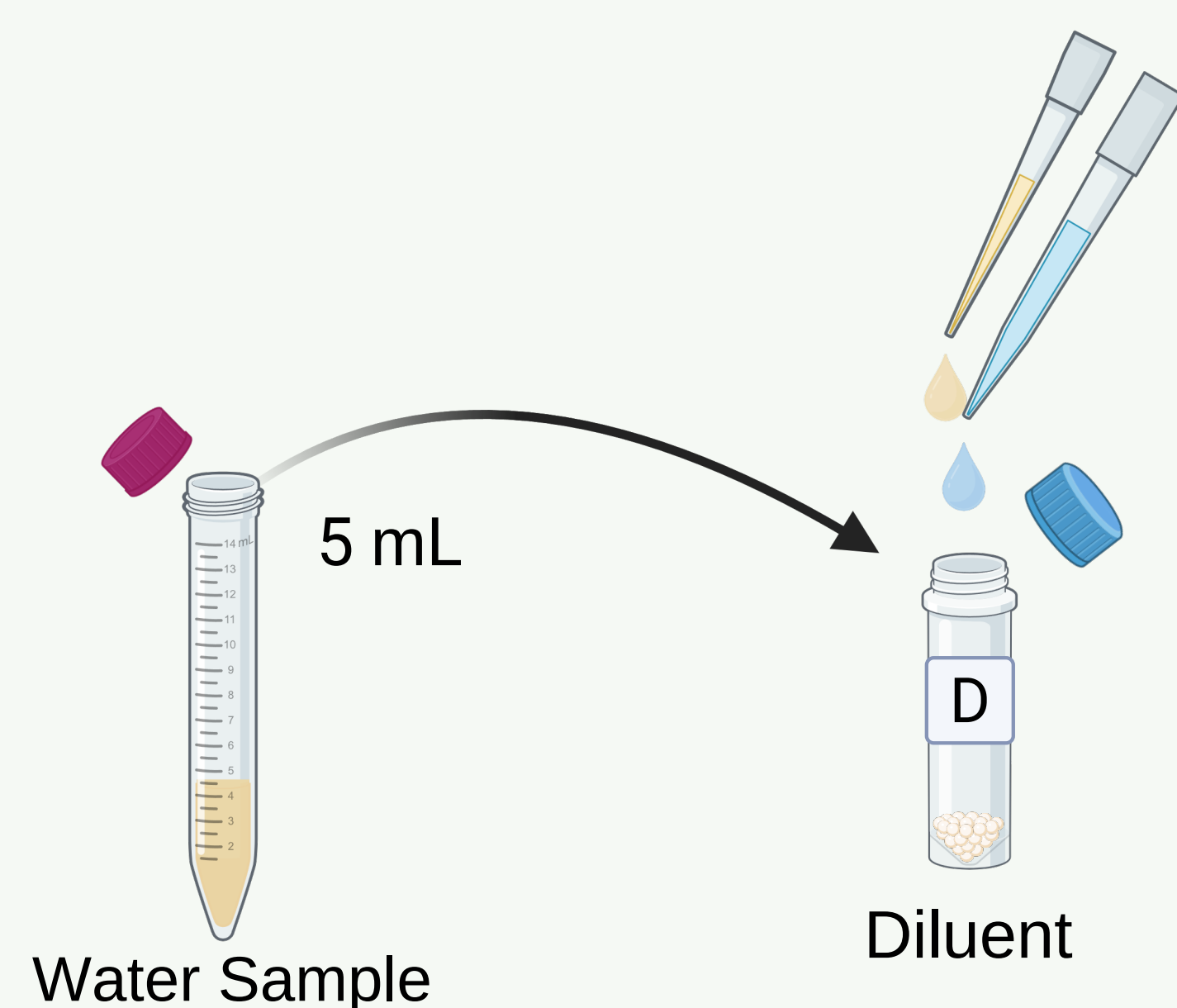

2. Add **140 µL** of **mixed sample** to a Conjugate tube. Pipette up and down to dissolve. Let sit uninterrupted at room temperature for 10 minutes.

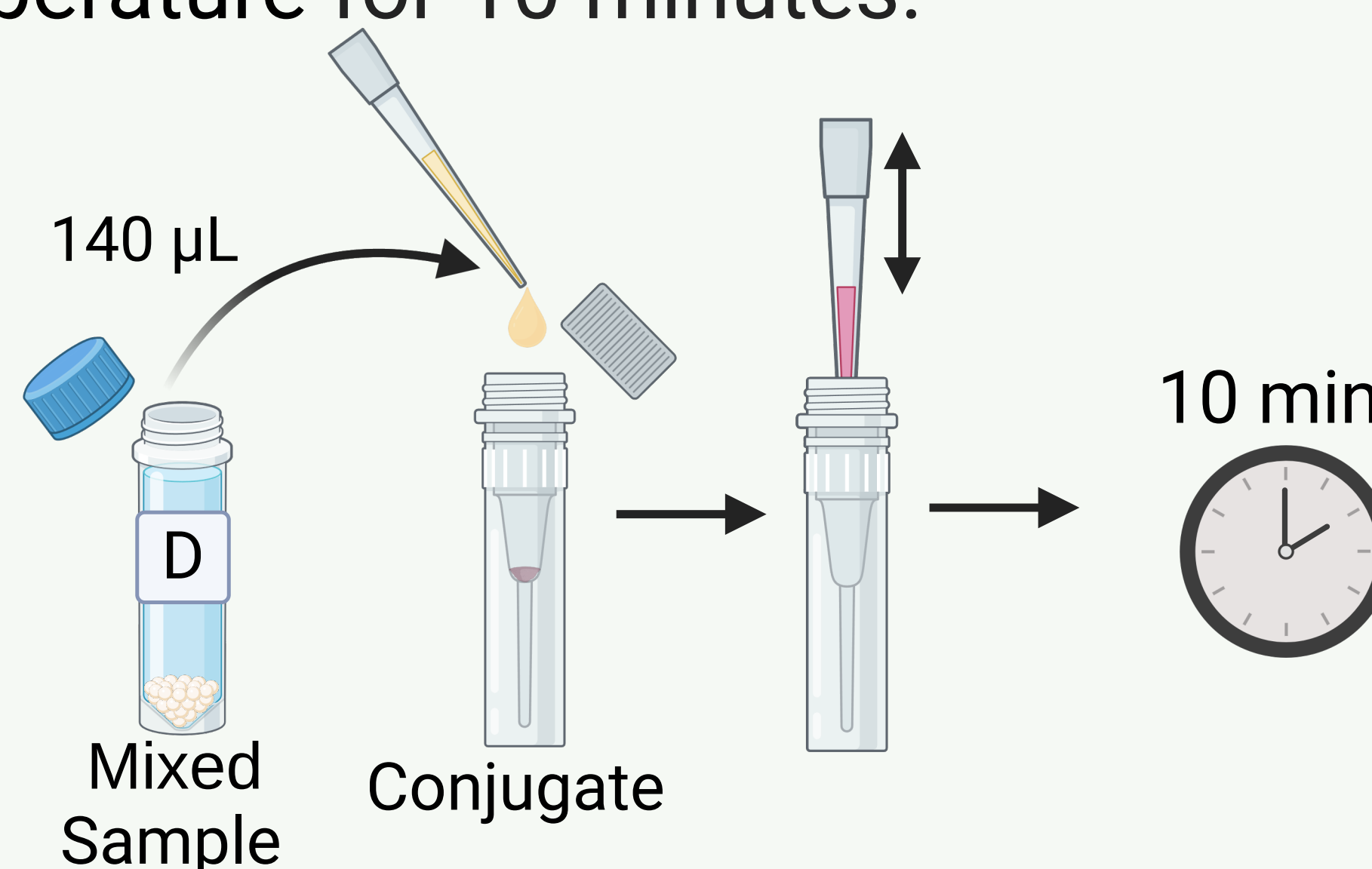


3. Add **all of the Conjugate mix solution** to the test cassette. Leave the cassette on a horizontal surface and set a timer for 15 minutes.

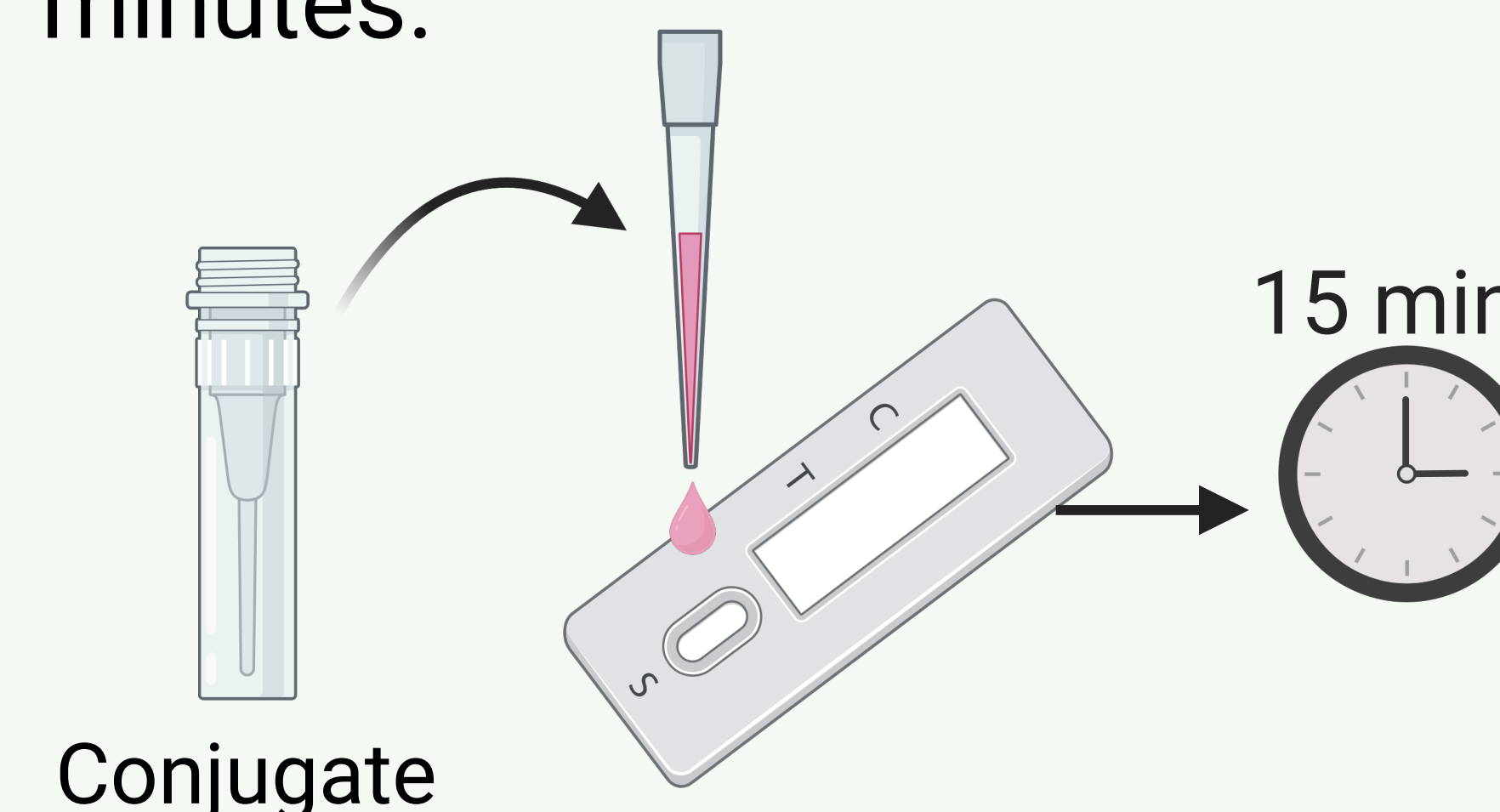


PART 2: RESULTS

Insert the cassette into the reader slot

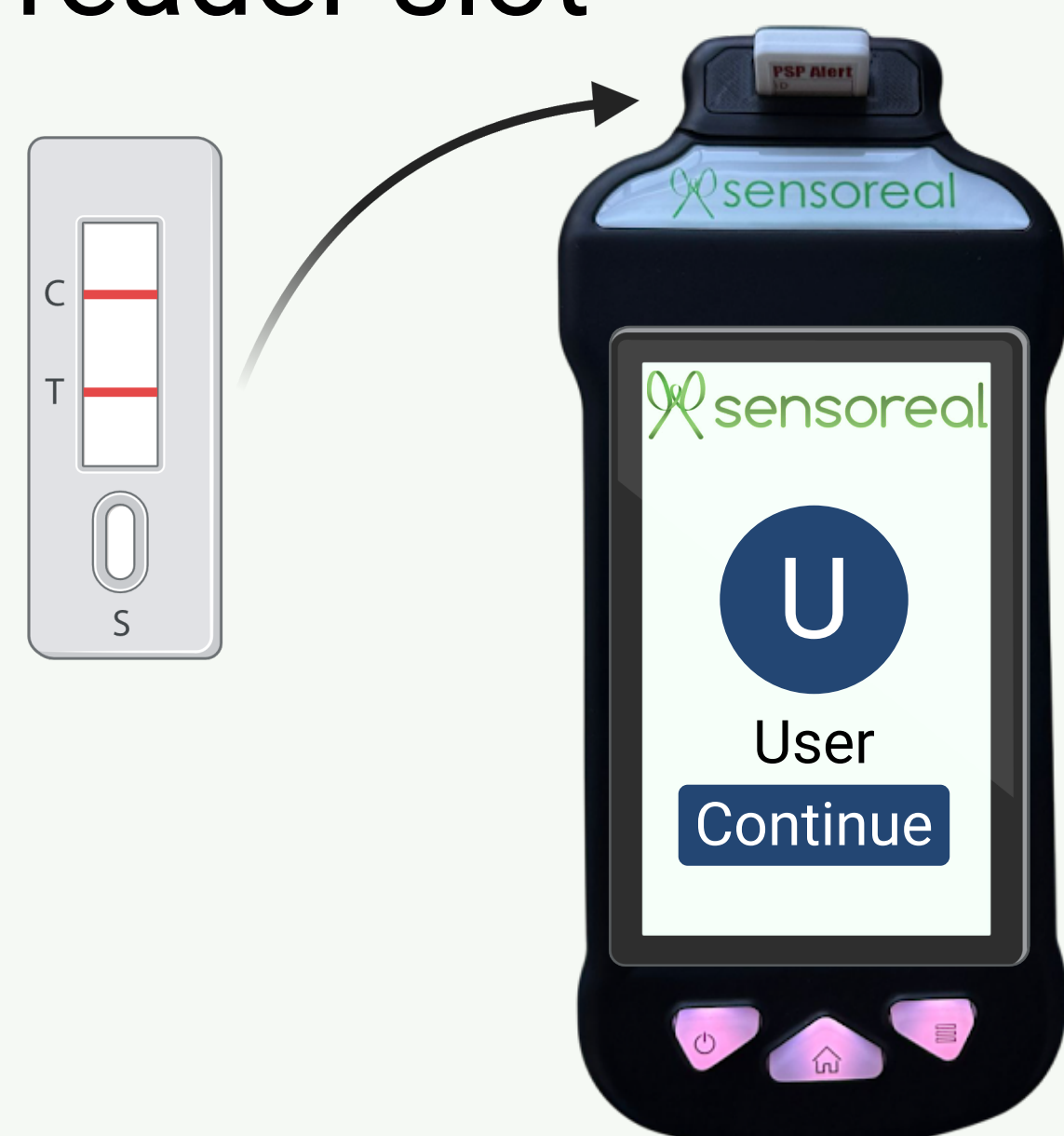

Follow the steps on the screen

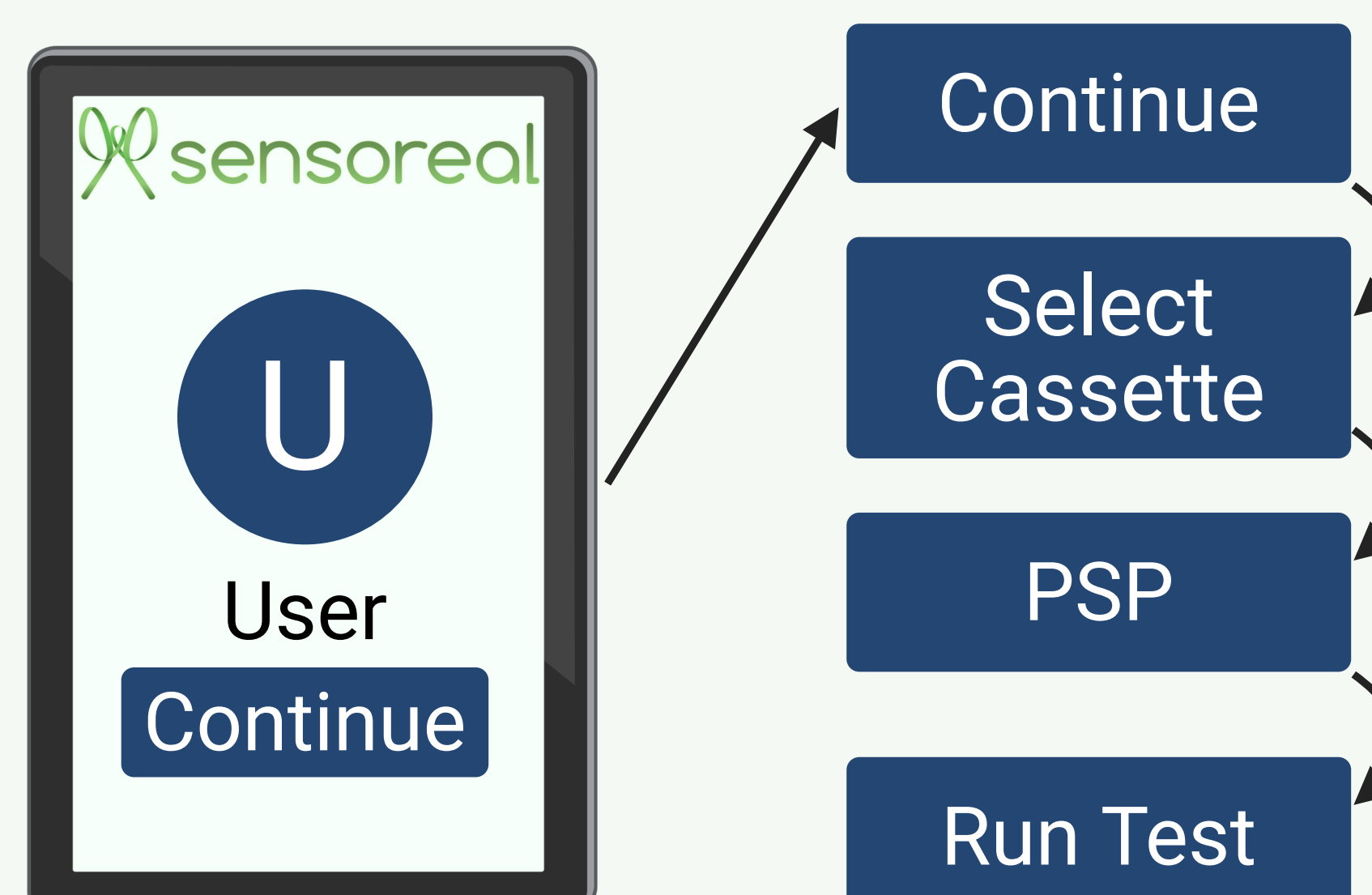

Test results will be displayed after a few seconds

